

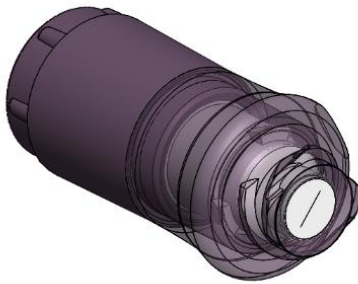
**DESCRIPTION**

FlowArt® BM01TY is a one-way connector with an integrated silicone seal valve.

**FEATURES**

- ✓ Microbial barrier swabable silicon stem with a split septum, flat seal to allow enhanced disinfection and drying time.
- ✓ Clear transparent housing to allow clear view into the fluid pathway for enhanced flushing.
- ✓ Allows flow only towards the patient.
- ✓ Compatible with TPN solutions, blood and oncology drugs.
- ✓ Suitable for high pressure infusions (up to 325 psi).
- ✓ Allows multiple accesses.
- ✓ ISO 80369-7 compatible.
- ✓ Photoprotective (-PH) versions are in amber color.
- ✓ Latex-free, PVC-free, DEHP-free

Part	Material
Valve Body	Copolyester
Swabable septum	Silicone
Check Valve Body	Silicone

**DIAGRAM****PRODUCT DATA**

Characteristic	Specification
Product Code	BM01TY, BM01TY-PH,
Product Name	Flowart® One-Directional Needle-Free Valve
Priming Volume (mL)	0.13
Displacement	0.00568 mL (semi-neutral)
Flow Rate (mL)	114
Flushing Volume (mL)	2
Length (cm)	3

**LABELLING AND PACKAGING**

Units are double bagged with polyethylene bags, placed in corrugated shipper, labelled on its bag and packaging container. Each container is labelled with the following traceability information:

- Quantity,
- Product description,
- Product date
- Lot number

TDS.108 \_Rev03\_EZ\_26.04.2021

ASSET® MEDICAL DESIGN INC is an ISO 13485 registered company, Certificate No.: M 8072.

FlowArt is a registered trademark of Asset Medikal.

©2021 Asset Medikal. All rights reserved.

Asset Medikal Tasarım Sanayi ve Ticaret A.Ş.

İkitelli OSB Mahallesi  
17. Cadde, No: 17  
Başakşehir, İstanbul - Turkey

Phone : +90 212 494 27 27

Fax : +90 212 494 34 37

Taxation Authority / No : İKİTELLİ V.D. / 0910501228

Trade Register No : 695458

Central Registration System No : 0091050122800019

info@assetmedikal.com | www.assetmedikal.com

**FlowArt®**

**asset®**  
Design the future