

## NEUTRART® NEEDLE-FREE VALVE, NEUTRAL

## DESCRIPTION


NeutArt® BM01N is a connector with an integrated silicone seal valve.

## FEATURES

- ☑ Microbial barrier swabable silicon stem with flat seal to allow enhanced disinfection and drying time.
- ☑ Clear transparent housing to allow clear view into the fluid pathway for enhanced flushing.
- ☑ Compatible with TPN solutions, blood and oncology drugs.
- ☑ Allows multiple accesses.
- ☑ ISO 80369-7 compatible.
- ☑ MRI and CT suitable.
- ☑ Photoprotective (-PH) versions are in amber color.
- ☑ Latex-free, PVC-free, DEHP-free

Part	Material
Valve Body	Copolyester
Swabable septum	Silicone

## DIAGRAM

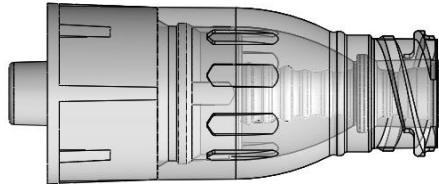

## PRODUCT DATA

Characteristic	Specification
Product Code	BM01N, BM01N-PH
Product Name	Needle-free Valve, Neutral
Priming Volume (mL)	0.04
Displacement (mL)	<0.01 (neutral)
Flow Rate (mL/min @ 1m water head pressure)	160
Flushing Volume (mL)	0.2
Length (cm)	2.5

## LABELLING AND PACKAGING

Each unit is labelled on its blister and packaging container.

Each container is labelled with the following traceability information:

- Quantity,
- Product description,
- Product date,
- Lot number

TDS.113 \_Rev02\_26.04.2021

ASSET® MEDICAL DESIGN INC is an ISO 13485 registered company, Certificate No.: M 8072.

NeutrArt is a registered trademark of Asset Medikal.

Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions.

©2021 Asset Medikal. All rights reserved.

## Asset Medikal Tasarım Sanayi ve Ticaret A.Ş.

İkitelli OSB Mahallesi  
17. Cadde, No: 17  
Basaksehir, Istanbul - Turkey

Phone : +90 212 494 27 27

Fax : +90 212 494 34 37

Taxation Authority / No : İKİTELLİ V.D. / 0910501228

Trade Register No : 695458

Central Registration System No : 0091050122800019

info@assetmedikal.com | www.assetmedikal.com